



***Innovating Health & Wellness  
Within Your World***

**ERGONOMIC ADVANTAGE**  
ULTRASOUND TABLES  
& ACCESSORIES

- Echocardiography
- Vascular
- Multi-Specialty
- Women's Imaging
- General Ultrasound



# ERGONOMIC ADVANTAGE ULTRASOUND TABLES

The new **EA** Line of tables incorporates all of the ergonomic features and comfort our tables have had for years with new powered features that speed up exam times and make positioning patients easier. From powered cardio panels and access panels to multiple motions occurring at the same time, these tables live up to their name.

Inspired by the leading ergonomic experts in sonography




**SOUND | ERGONOMICS**

- Decrease sonographer injury
- Minimize exam times



## TABLE OF CONTENTS

Echocardiography EA	4,5
Vascular with Fowler EA	6, 7
Multi-Specialty EA	8, 9
Womens Imaging EA	10, 11
General 3 Section EA	12, 13
General 2 Section EA	14, 15
Features and Options	16
Accessories	17

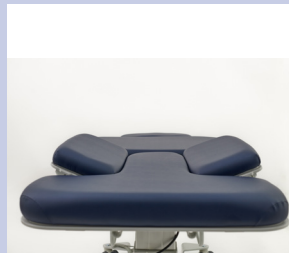


# ERGONOMIC ADVANTAGE ULTRASOUND TABLES

## FOCUSED ON ERGONOMICS

Designed to reduce Repetitive Stress Injuries by minimizing reaching and the energy required to position patients. Our new designs support patients better, improve access, and make it faster to position patients.



Large Access panel supports more client anatomy while minimizing reaching for sonographers.

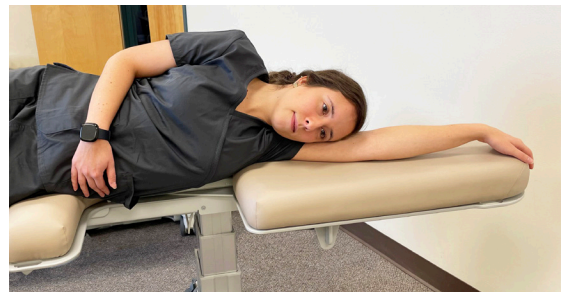

Powered Echocardiography panel slides open quickly to the desired size, simultaneously balancing access with support for the patient's torso and arm.

---

## DESIGNED TO LAST

Made in America, all of our tables feature our unique small profile base made of heavy duty, powder coated steel and twin lifting towers for the ultimate in strength and rigidity.



Replaceable padded sections that are easily removed.




Powder coated, heavy steel construction with twin lifting towers that accommodate patients up to 550lbs



## Ergonomic Advantage Line ECHOCARDIOGRAPHY EA ULTRASOUND TABLE



Shown in Clay fabric

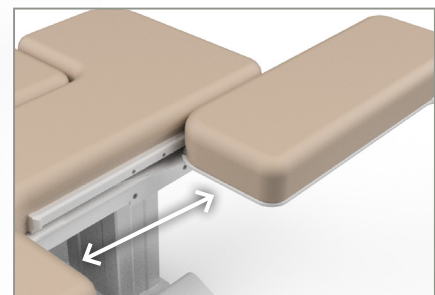

Sliding Cardio Panel  
opens as wide as 15" (38cm)

Inspired by the leading ergonomic experts in sonography



- Decrease sonographer injury
- Minimize exam times

> **View Video**




### Recommended for:

- Cardiac exams
- Vascular exams
- All abdominal & superficial structures exams
- Musculoskeletal exams
- Biopsies

# ECHOCARDIOGRAPHY EA

## Fully Electric Table

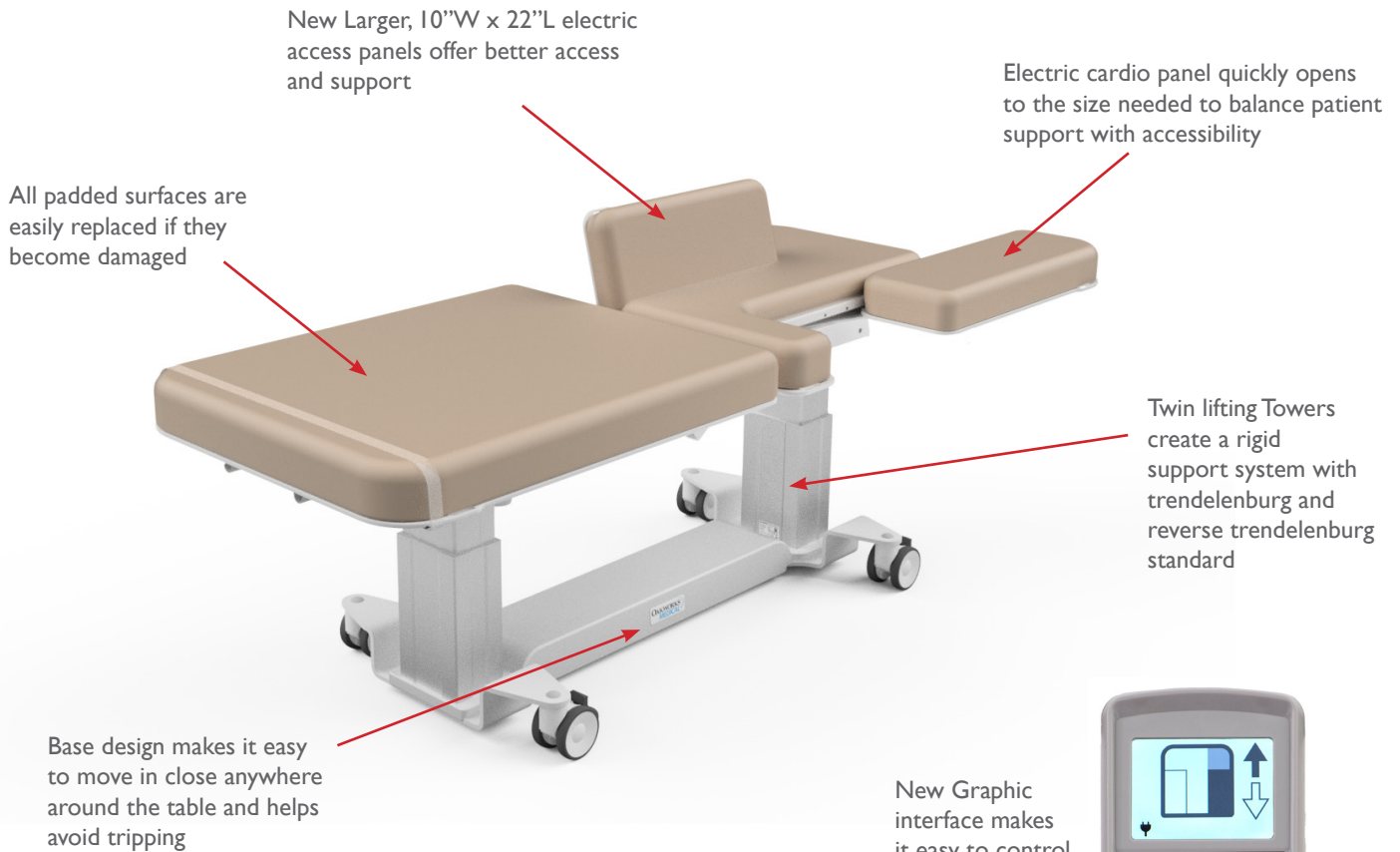

## ECHOCARDIOGRAPHY TABLE EA

#84773-T\*

FEATURES	
2 Section Top	30" (76cm) wide x 72" (183cm) long
Backrest	34" (86cm) long
Seat	38" (97cm) long
Motion 1	Electric Height Range 19" - 33" (48-83cm)
Motion 2	Electric Trendelenburg/Reverse Tilt $\pm 18^\circ$
Motion 3	Electric Fowler Backrest 0-75°
Motion 4	Electric Sliding Cardio Panel opens 15" (38cm)
Motion 5	Electric Access Panel opens 10" wide (25cm) x 0" - 22" long (56cm)
Hand Control	All electric motions
Weight Capacity	550 lb (250kg)
Product Weight	359 lb (163kg)
Padding	3" (8cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Included Features	4" (10cm) locking casters, Powered access back support, electric sliding cardio panel, Paper roll holder
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, CE marked, cETLus listed, CB Certificate

OPTIONS <i>must be ordered with table.</i>	
Central Base Locking System	Powered, lock and unlock using hand control
Foot Control	Height, Trendelenburg, Fowler
Side Rails (set of 2)	25.5" (65cm) long x 6" (15cm) tall (from surface of table)
Retractable footrests	Stirrups - adjust in length and width
Procedure tray	12.25" (6cm) x 10.25" (26cm) and holds 2 quarts (1.9 liters)
Battery Backup	Capable of running during power failure
Non-USA electric system	please contact us for order details

ACCESSORIES <i>order any time.</i>		#
Detachable Arm Board	Use with T-rail adapter 5" (13cm) wide x 20.5" (52cm) long	84943-T*
T-rail Adapter 8" (20cm) long	Attaches to the table for use with USA format T-rail attachable accessories.	83139
Table Extender	24" (61cm) wide x 12" (30cm) long	63487-T*
QuickLock™ Head Rest & AeroCel Pad	Comfortable, adjustable prone or supine support	PKG7554-T*

\*Fabric color code must be specified.



# Ergonomic Advantage Line VASCULAR WITH FOWLER EA ULTRASOUND TABLE



“Finally we have a table designed for vascular ultrasound!”

Patricia A Poe, BA RVT FSVU

Director Noninvasive Vascular Laboratory Services, UC Health



**SOUND | ERGONOMICS**

Inspired by the leading ergonomic experts in sonography



Shown in Ruby fabric

> **View Video**




## Recommended for:

- Venous reflux studies
- All other vascular exams
- Abdominal & superficial structures exams
- Musculoskeletal exams
- Biopsies

# VASCULAR with Fowler EA

## Fully Electric Table

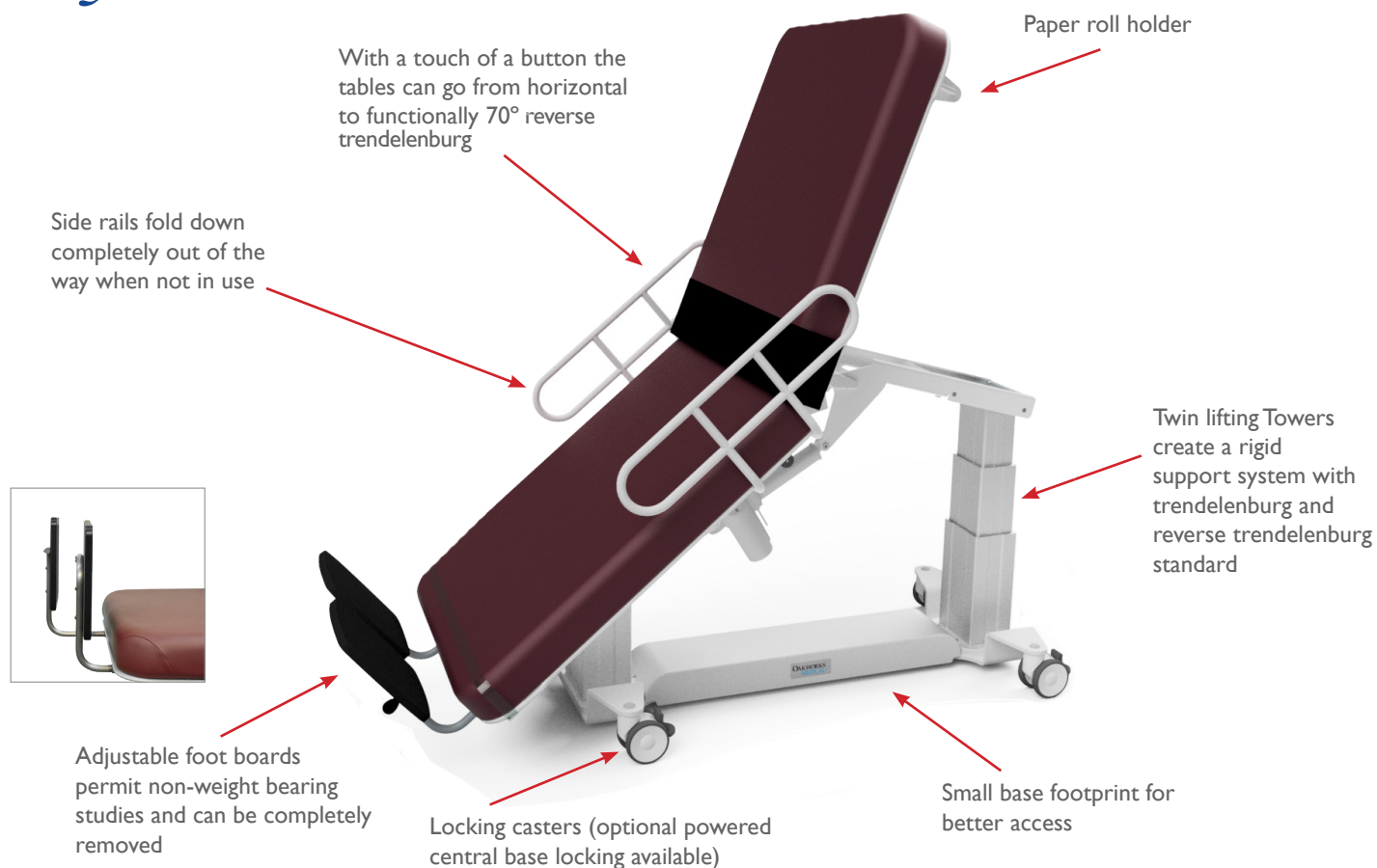

## VASCULAR WITH FOWLER TABLE EA

#84792- 26" W #83917 - 30" W

FEATURES	
2 Section Top	26" (66cm) or 30" (76cm) wide x 72" (183cm) long
Backrest	34" (86cm) long
Seat	38" (97cm) long
Motion 1	Electric Height Range 19" - 33" (48-83cm)
Motion 2	Electric Trendelenburg/Reverse Tilt ±18°
Motion 3	Electric Additional Reverse Trendelenburg 0-25° (40° total)
Motion 4	Electric 0-75° Backrest (at 40° max tilt, 70° functional)
Hand Control	All powered motions
Weight Capacity	550 lb (250kg)
Product Weight	359 lb (163kg)
Padding	3" (8cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Included Features	4" (10cm) locking casters, Integrated foot board (split, drop away), Paper roll holder, Patient safety strap, Side rails (2) 25" (64cm) long x 6" (15cm) tall (from surface of table)
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, CE marked, cETLus listed, CB Certificate

OPTIONS must be ordered with table.	
Cardio Panel and Powered Access/Back Support	Electric Sliding Cardio Panel opens 15" (38cm) - left side Electric Access Panel opens 10" wide (25cm) x 22" long (56cm)- right side
Powered Central Base Locking	Lock and unlock using hand control
Foot Control	Height, Trendelenburg, Backrest
Battery Backup	Capable of running during power failure
Non-USA Electric System	Please contact us for order details

ACCESSORIES order any time.		#
Detachable Arm Board	Use with T-rail adapter 5" (13cm) wide x 20.5" (52cm) long	84943-T*
T-rail Adapter 8" (20cm) long	Attaches to the table for use with USA format T-rail attachable accessories.	83139
Table Extender	24" (61cm) wide x 12" (30cm) long	63487-T*
QuickLock™ Head Rest & AeroCel Pad	Comfortable, adjustable prone or supine support	PKG7554-T*

\*Fabric color code must be specified.



## Ergonomic Advantage Line MULTI-SPECIALTY EA ULTRASOUND TABLE



Shown in  
Sapphire fabric



Shown with stirrups

Inspired by the leading ergonomic experts in sonography



**SOUND | ERGONOMICS**

- Decrease sonographer injury
- Minimize exam times

> **View Video**



### *Recommended for:*

- All abdominal & superficial structures exams (*thyroid, male pelvis*)
- Musculoskeletal exams
- OB/GYN exams
- Breast exams
- Vascular exams
- Biopsies



# MULTI-SPECIALTY EA

## Fully Electric Table

New Powered Leg section makes positioning patients' legs for exams stress free

All padded surfaces are easily replaced if they become damaged

New Side Rails are recessed even farther under the top when not in use (optional)

New Larger, 10"W x 22"L access panels are powered and offer better access and support

Paper roll holder

Twin lifting Towers create a rigid support system with trendelenburg and reverse trendelenburg standard

New Electric Central Base Locking system makes securely positioning the table fast and easy (optional)

Base design makes it easy to move in close anywhere around the table on a chair

New Graphic interface makes it easy to control all the powered functions.



### MULTI-SPECIALTY TABLE EA

#84797 30" W

FEATURES	
3 Section Top	30" (76cm) wide x 72" (183cm) long
Backrest	34" (86cm) long
Seat	26" (66cm) long
Leg Rest	11" (28cm) long x 19" (48cm) wide
Motion 1	Electric Height Range 19" - 33" (48-83cm)
Motion 2	Electric Trendelenburg/Reverse Tilt ±18°
Motion 3	Electric Fowler Backrest 0-75°
Motion 4	Electric Leg Rest 0-90°
Access/Back Support	Electric 10" wide (25cm) X 22" long (56cm) 0-75° Access/Back support for patient positioning (l & r side)
Hand Control	Graphic display with presets
Weight Capacity	550 lb (250kg)
Product Weight	359 lb (163kg)
Padding	3" (8cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Included Features	Stirrups - adjust in length and width, Powered access back support (left & right) 0-75°, 4" (10cm) locking casters, Paper roll holder
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, CE marked, cETLus listed, CB Certificate

OPTIONS <i>must be ordered with table.</i>	
Foot Control	Height, Trendelenburg/Reverse, Fowler, Leg Section
Procedure Tray	12.25" (6cm) x 10.25" (26cm) and holds 2 quarts (1.9 liters)
Side Rails (set of 2)	25.5" (65cm) long x 6" (15cm) tall (from surface of table)
Powered Central Base Locking	Lock and unlock using hand control
Battery Back Up	Capable of running during power failure
Non-USA electric system	Please contact us for order details

ACCESSORIES <i>order any time.</i>		#
Detachable Arm Board	Use with T-rail adapter (sold separately) 5" (13cm) wide x 20.5" (52cm) long	84943
QuickLock™ Head Rest & AeroCel Pad	Comfortable, adjustable prone or supine support	PKG7554-T*
Table Extender	24" (61cm) wide x 12" (30cm) long	63487-T*
T-rail Adapter 8" (20cm) long	Attaches to the backrest for use with USA format T-rail attachable accessories.	83139

\*Fabric code must be specified when ordering.



## Ergonomic Advantage Line WOMENS IMAGING EA ULTRASOUND TABLE



Shown in Forest fabric



Inspired by the leading ergonomic experts in sonography



**SOUND | ERGONOMICS**

- Decrease sonographer injury
- Minimize exam times


> **View Video**



Recommended for:

- OB/GYN exams
- All abdominal & superficial structures exams
- Musculoskeletal exams
- Vascular exams
- Biopsies

# WOMENS IMAGING EA

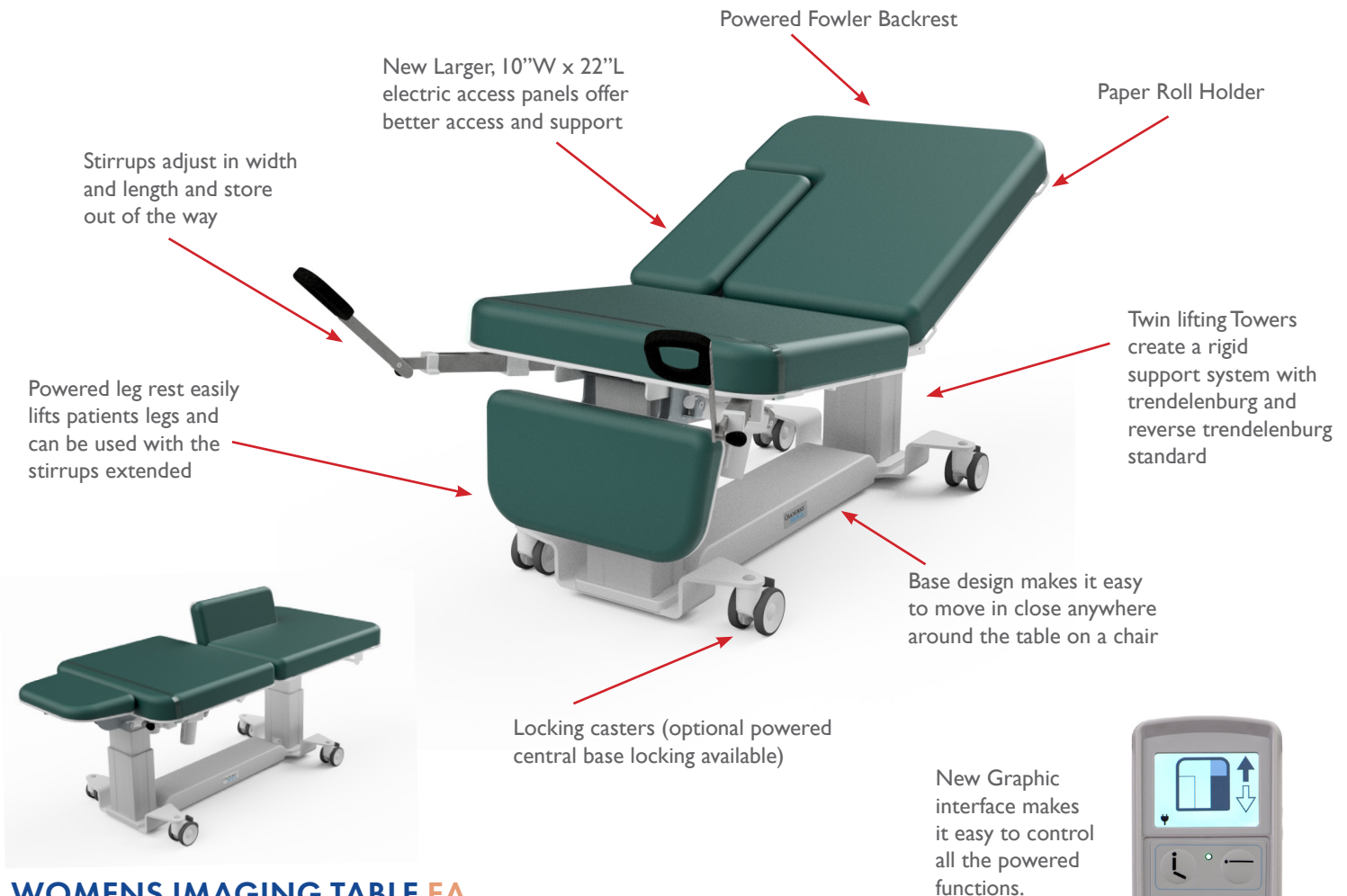

## WOMENS IMAGING TABLE EA

#84472-T\* 30" W

FEATURES	
3 Section Top	30" (76cm) wide x 72" (183cm) long
Backrest	34" (86cm) long
Seat	26" (61cm) long
Leg Rest	11" (28cm) long x 19" (48cm) wide
Motion 1	Electric Height Range 19" - 33" (48-83cm)
Motion 2	Electric Trendelenburg/Reverse Tilt $\pm 18^\circ$
Motion 3	Electric Fowler Backrest 0-75°
Motion 4	Electric Leg Rest 0-90°
Access/Back Support	Electric 10" wide (25cm) X 22" long (56cm) 0-75° Access/Back support for patient positioning (right side)
Hand Control	All powered motions
Weight Capacity	550 lb (250kg)
Product Weight	359 lb (163kg)
Padding	3" (8cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Included Features	4" (10cm) locking casters, Paper roll holder, Stirrups (Retractable foot rests)
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, CE marked, cETLus listed, CB Certificate

OPTIONS must be ordered with table.	
Foot Control	Height, Trendelenburg, Fowler
Powered Central Base Locking	Lock and unlock using hand control
Procedure Tray 30" (76 cm)	12.25" (6cm) x 10.25" (26cm) and holds 2 quarts (1.9 liters)
Side Rails (set of 2)	25.5" (65cm) long x 6" (15cm) tall (from surface of table)
Battery Backup	Capable of running during power failure
Non-USA Electric System	Please contact us for order details

ACCESSORIES order any time.		#
Detachable Arm Board	Use with T-rail adapter 5" (13cm) wide x 20.5" (52cm) long	84943
QuickLock™ Head Rest (black) & AeroCel Pad	Comfortable, adjustable prone or supine support	PKG7554*
Table Extender	24" (61cm) wide x 12" (30cm) long	63487-T*
T-rail Adapter 8" (20cm) long	Attaches to the backrest for use with USA format T-rail attachable accessories.	83139

\*Fabric color code must be specified.



## Ergonomic Advantage Line GENERAL 3 SECTION EA ULTRASOUND TABLE



Shown in Pewter fabric




Inspired by the leading ergonomic experts in sonography



- Decrease sonographer injury
- Minimize exam times

> **View Video**

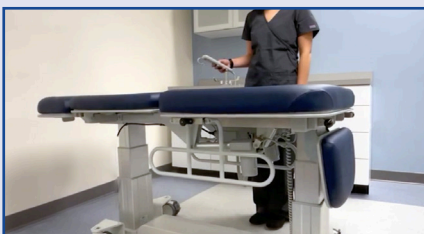

Visit  
[www.oakworksmed.com](http://www.oakworksmed.com)  
for a quote!

# GENERAL 3 SECTION EA



## GENERAL 3SECTION TABLE EA


#84783-T\* 26" W    #83892-T\* 30" W

FEATURES	
3 Section Top	26" (66cm) or 30" (76cm) wide x 72" (183cm) long
Backrest	34" (86cm) long
Seat	26" (66cm) long
Leg Rest	19" (48cm) wide x 11" (28cm) long
Motion 1	Electric Height Range 19" - 33" (48-84cm)
Motion 2	Electric Trendelenburg/Reverse Tilt ±18°
Motion 3	Electric Fowler Backrest 0-75°
Motion 4	Electric Leg Rest 0-90°
Controls	Multi-function hand control
Weight Capacity	550 lb (250kg)
Product Weight	359 lb (163kg)
Padding	3" (8cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Included Features	4" (10cm) locking casters, Paper roll holder
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, CE marked, cETLus listed, CB Certificate

OPTIONS must be ordered with table.	
Foot Control	Height, Trendelenburg, Fowler
Powered Central Base Locking	Lock and unlock using hand control
Procedure Tray	12.25" (6cm) x 10.25" (26cm) and holds 2 quarts (1.9 liters)
Retractable footrests	Stirrups - adjust in length and width
Side Rails (set of 2)	25.5" (65cm) long x 6" (15cm) tall (from surface of table)
Battery Backup	Capable of running during power failure
Non-USA Electric System	Please contact us for order details

ACCESSORIES order any time.		#
Detachable Arm Board	Use with T-rail adapter 5" (13cm) wide x 20.5" (52cm) long	84943
QuickLock™ Head Rest (black) & AeroCel Pad	Comfortable, adjustable prone or supine support	PKG7554*
Table Extender	24" (61cm) wide x 12" (30cm) long	63487-T*
T-rail Adapter 8" (20cm) long	Attaches to the backrest for use with USA format T-rail attachable accessories.	83139

\*Fabric color code must be specified.

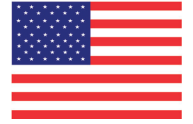

## Ergonomic Advantage Line GENERAL 2 SECTION EA ULTRASOUND TABLE



Shown in Ruby fabric

Inspired by the leading ergonomic experts in sonography



**SOUND | ERGONOMICS**


- Decrease sonographer injury
- Minimize exam times

> **View Video**



Visit  
[www.oakworksmed.com](http://www.oakworksmed.com)  
for a quote!

# GENERAL 2 SECTION EA

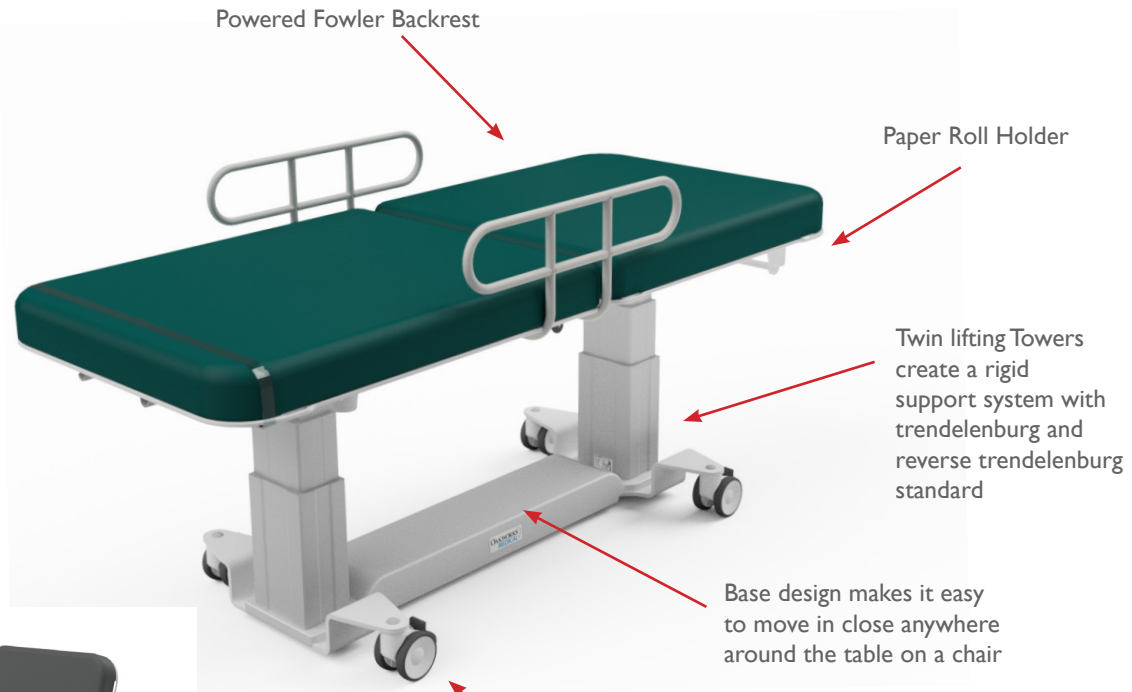

Locking casters' double wheel design makes positioning the table easy

New Graphic interface makes it easy to control all the powered functions.



## GENERAL 2 SECTION TABLE EA

#84771-T\* 26" W #83901-T\* 30" W

FEATURES	
2 Section Top	26" (66cm) or 30" (76cm) wide x 72" (183cm) long
Backrest	34" (86cm) long
Seat	38" (97cm) long
Motion 1	Electric Height Range 19" - 33" (48-84cm)
Motion 2	Electric Trendelenburg/Reverse Tilt $\pm 18^\circ$
Motion 3	Electric Fowler Backrest 0-75°
Controls	Multi-function hand control
Weight Capacity	550 lb (250kg)
Product Weight	359 lb (163kg)
Padding	3" (8cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Included Features	4" (10cm) locking casters, Paper roll holder
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, CE marked, cETLus listed, CB Certificate

OPTIONS must be ordered with table.	
Foot Control	Height, Trendelenburg, Fowler
Powered Central Base Locking	Lock and unlock using hand control
Procedure Tray 30" (76 cm)	12.25" (6cm) x 10.25" (26cm) and holds 2 quarts (1.9 liters)
Side Rails (set of 2)	25.5" (65cm) long x 6" (15cm) tall (from surface of table)
Retractable footrests	Stirrups - adjust in length and width
Battery Backup	Capable of running during power failure
Non-USA Electric System	Please contact us for order details

ACCESSORIES order any time.		#
Detachable Arm Board	Use with T-rail adapter 5" (13cm) wide x 20.5" (52cm) long	84943
QuickLock™ Head Rest (black) & AeroCel Pad	Comfortable, adjustable prone or supine support	PKG7554*
Table Extender	24" (61cm) wide x 12" (30cm) long	63487-T*
T-rail Adapter 8" (20cm) long	Attaches to the backrest for use with USA format T-rail attachable accessories.	83139

\*Fabric color code must be specified.

# ERGONOMIC ADVANTAGE LINE OPTIONS

## SIDE RAILS

Our new design recesses farther under the top when not in use, improving ergonomics beyond our previous model.



## POWERED CENTRAL BASE LOCKING

The new system allows you to easily lock or unlock the base by simply pushing a button on the hand control.



## FOOT CONTROL

The foot control allows you to quickly position the patient by controlling multiple functions without using your hands.




## BATTERY BACK UP

A battery back up system capable of running for several hours if the power fails.



## STIRRUPS

Easily adjusted stirrups that store under the top can be combined with our optional, retractable, removable procedure tray.

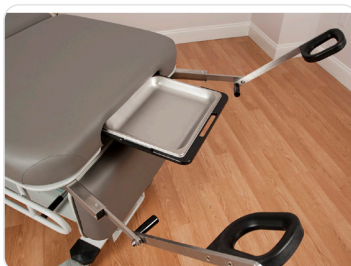

## PROCEDURE TRAY

Retractable, removable stainless steel. 12.25" (6cm) x 10.25" (26cm) & holds 2 quarts (1.9 liters).





# ERGONOMIC ADVANTAGE LINE ACCESSORIES



## SONOGRAPHER CHAIRS

Sonography chairs provide adjustments for seat height, back/arm height, and backward tilt tension, along with an integrated seat slider. These adjustments, along with the ergonomically designed support surfaces, allow sonographers to find the least stressful position for a wide variety of scans.

## T-RAIL ADAPTER

Multiple locations allow you to choose the best spot for this easily installed & removed accessory capable of handling up to 150 lb (68kg). Can be used with any American standard accessory. 8" (20cm) L



## TABLE EXTENDER

Lengthens the table by 12" (31cm) to 84" (213cm) allowing for access to the patient's neck for carotid artery exams and thyroid exams. 24" (61cm) W x 12" (31cm) L



## DETACHABLE ARM BOARD

Attaches to the T-rail Adapter and rotates over 140 degrees. 5" (13cm) wide x 20.5" (52cm) long (Requires T-Rail Adapter)



## QUICKLOCK™ HEAD REST & AERO · CEL PAD


The platform allows precise cervical flexion and height positioning over a wide range for both prone and supine treatments.



# Notes

# Notes

# TERRATOUGH™ FABRIC

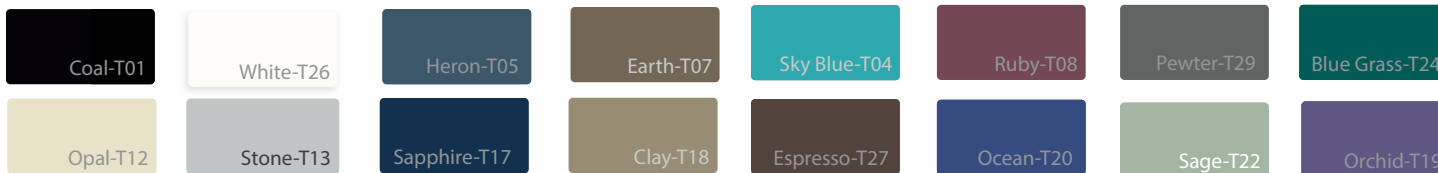

*We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible. However, due to the limitations of printers, scanners, computer monitors, and other equipment, the colors may not be accurate.*



Oakwork's is independently certificated to ISO 13485 Quality Management System and regulatory requirement for medical devices ensuring consistent, high quality products. Our products are made in our state-of-the-art 91,000 square foot factory in New Freedom, Pennsylvania.



All of our tables are FDA Listed.

Please visit our website, [www.oakworksmed.com](http://www.oakworksmed.com), to view our complete line of tables and accessories for C-arm Imaging, Pain Management, Ultrasound, Urology, Lithotripsy, Exam & Treatment.



Call For More Information

Oakworks Medical 717.759.3146 [www.oakworksmed.com](http://www.oakworksmed.com)