

OAKWORKS<sup>®</sup>  
**MEDICAL**  
*Position for Success*

# CFPMB301

*C-Arm Imaging Table*

- **Safety Tested and Approved**  
FDA listed, CE marked,  
cETLus approved, CB certificate
- **Industry Leading Warranty**  
3 years parts, 2 years labor
- **Superior Height Range**  
super low 26" and high 44"  
*Great for patient transfers!*

*Ask us about  
our 5 Day  
Quick Ship  
Lead Time  
Program!  
does not include  
transit time*



**750 LB.  
BARIATRIC  
PATIENT WEIGHT**



ISO 13485:2003  
FM 584901

made in the USA  
with US & imported parts



**Contact Your Oakworks Factory Rep**

FEATURES	
Carbon Fiber Top 52.5" (133cm) Metal free imaging space	Rectangular 24" (61cm) wide x 80" (203cm) long
Motion 1	Electric Height Range 26" - 44" (66-112cm)
Motion 2	Electric Trendelenburg tilt +15°/-11°
Motion 3	Electric Longitudinal travel 4" (10cm)
Controls	Multi-function hand control
Weight Capacity	750 lb (340kg)
Product Weight	475 lb (216kg)
Included Features	1" (3cm) Radiolucent Comfort Foam™ table top pad, 5" (13cm) locking casters, patient safety straps (2)
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, CE marked, cETLus approved, CB Certificate

OPTIONS <i>must be ordered with table.</i>		#
2" (5cm) Table Top Pad Upgrade	upgrade from the standard 1" (3cm) pad	73190-T01
Foot Control	multi-functional	75076
Non-USA electric system	please contact us for order details	73134

ACCESSORIES <i>order any time.</i>		#
Carbon Fiber Arm Board (6" wide)	180° Articulation	42585-T01
Carbon Fiber Arm Board (7.5" wide)	180° Articulation	70337-T01
Fluoro Extender	Polycarbonate	33975
Spine Positioning System II	8 piece system	70100-T01
T-rail Clamp	6" (15cm) long rail	63336

## Table Movements



26" - 44"  
HEIGHT RANGE



4"  
LONGITUDINAL TRAVEL




+15°/-11°  
TRENDLENBURG TILT

## ACCESSORIES

### CARBON FIBER ARM BOARD

- 180° range of motion
- Requires no attachment hardware, slides under table top pad.

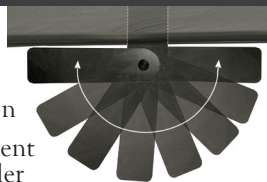

#42585-T01 - 6" wide  
#70337-T01 - 7.5" wide

### SPINE POSITIONING SYSTEM II

- Optimal patient positioning.
- Improved imaging.



#70100-T01

### 6" ADJUSTABLE T-RAIL CLAMP

- Accommodates a variety of T-Rail attachable accessories.



#63336

**OAKWORKS<sup>®</sup>**  
**MEDICAL**  
*Position for Success*

View our full line of  
Medical products



Scan this code with  
your smart phone

#72703



717.235.6807 [www.oakworksmed.com](http://www.oakworksmed.com) email: [customerservice@oakworks.com](mailto:customerservice@oakworks.com)