

# GENERAL WITH FOWLER Ultrasound Table



Start with a General Ultrasound Table, then simply add the features you want based on the type of ultrasound specialty you perform.

- 2 section top table with many optional features to customize your needs
- Dynamic 22" to 38" (56-97cm) adjustable height range for wheelchair transfer & ergonomic working heights
- 550 lb (250kg) bariatric rated lift & load capacity
- Twin tower design for added strength & rigidity



shown in Ruby fabric



Trendelenburg/Reverse



22" - 38" Height Range

Inspired by the leading ergonomic experts in sonography



**SOUND | ERGONOMICS**

- Decrease sonographer injury
- Features minimize exam times

Build-A-Bed™ is Trademarked by Sound Ergonomics.



FM584901

made in the USA  
with US & imported parts



**OAKWORKS®**  
**MEDICAL**  
*Position for Success*

# GENERAL WITH FOWLER Ultrasound Table

## Recommended for:

- All abdominal & superficial structures exams
- Musculoskeletal exams
- Vascular exams
- Biopsies



shown in Ruby fabric with optional Central Base Locking System and Side Rails

## GENERAL WITH FOWLER TABLE

#69222-T\* 27" W    #61988-T\* 30" W

FEATURES	
2 Section Top	27" (69cm) or 30" (76cm) wide x 72" (183cm) long
Backrest	34" (86cm) long
Seat	38" (97cm) long
Motion 1	Electric Height Range 22" - 38" (56-97cm)
Motion 2	Electric Trendelenburg/Reverse Tilt $\pm 15^\circ$
Motion 3	Electric Fowler Backrest 0-75°
Controls	Multi-function hand control
Weight Capacity	550 lb (250kg)
Product Weight	359 lb (163kg)
Padding	3" (8cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Included Features	4" (10cm) locking casters, Paper roll holder
Warranty	3 years parts, 2 years labor

\*Fabric color code must be specified.

OPTIONS must be ordered with table.		#
Accessory Mounting Plate	Fowler end	62123
Adjustable Access Back Support	Left side 0-85°	73133
Adjustable Access Back Support	Right side 0-85°	73132
Central Base Locking System	Not available with Patient Lift Compatible Base	61777
Foot Control	Height, Trendelenburg, Backrest	61901
Foot Control	Height, Trendelenburg, Backrest, Cardio Panel	65679
Non-USA electric system	please contact us for order details	73134
Patient Lift Compatible Base	Not available with Central Base Locking System	65893
Side Rails (set of 2)	23.5" (60cm) long x 6" (15cm) tall (from surface of table)	61781
Sliding Cardio Panel - Electric	30" (76cm) width table, left side opens 15" (38cm)	74080
Sliding Cardio Panel - Manual	30" (76cm) width table, left side opens 15" (38cm)	73130
ACCESSORIES order any time.		#
Detachable Arm Board	Use with T-rail adapter 5" (13cm) wide x 20.5" (52cm) long	62192
QuickLock™ Head Rest (black) & AeroCel Pad	Comfortable, adjustable prone or supine support (requires mounting plate)	PKG7554*
Table Extender	24" (61cm) wide x 12" (30cm) long (requires mounting plate)	63487

# ULTRASOUND FEATURES & OPTIONS

## CENTRAL BASE LOCKING SYSTEM #61777

Easily lock your table into place with our unique system that allows locking from both sides of the table.



OPTIONAL

## SIDE RAILS #61781

Increase stability and security for your patient. Fold down out of the way when not in use. 23.5" (60cm) long x 6" (15cm) tall (from surface of table)



OPTIONAL

## FULL FUNCTION FOOT CONTROLS

*#enquire for part number*

You can add a foot control that controls all functions.



*Multi-function foot control*

## PATIENT LIFT COMPATIBLE BASE #65893

This option allows the use of patient lifts and easy wheel chair transfers to minimize or eliminate caregiver injuries associated with patient transfers. The high clearance base permits all patient lifts to slide easily under the center of the table for stable transfers while not increasing the low height range of the tables that enable easy wheel chair transfers. Must be ordered with table.

- Includes larger, 5" (13cm) diameter Heavy Duty locking casters
- Provides 6.5" (17cm) clearance for patient lift
- Central Base Locking System not available with this option



OPTIONAL

## ADJUSTABLE ACCESS BACK SUPPORT

Left: #73133 Right: #73132

Adjustable with 5 positions 0-85°, allowing easy access for sonographer to the far side and back support for patients. 10" (25cm) wide x 15.5" (39cm) long



OPTIONAL

## SLIDING CARDIO PANEL *Electric: #74080 Manual: #73130*

The sliding cardio panel allows you to adjust the opening to the size you need, up to 15" (38cm). 10" (25cm) wide x 28" (71cm) long. Great for left or right handed scanning. Can also be used as an armrest. Available in manual or electric adjustment options. Available on left side of 30" (76cm) width tables.



OPTIONAL

## ACCESSORY MOUNTING PLATE

#62123

Add a face rest or table extender to your table.



OPTIONAL



# ULTRASOUND ACCESSORIES

## DETACHABLE ARM BOARD #62192-T\*

Attaches to the T-rail Adaptor and rotates over 140 degrees. 5" (13cm) wide x 20.5" (52cm) long (Requires T-Rail Adaptor)



## T-RAIL ADAPTOR #61902

Multiple locations allow you to choose the best spot for this easily installed & removed accessory capable of handling up to 150 lb (68kg). Can be used with any American standard accessory. 8" (20cm) long



## TABLE EXTENDER #63487-T\*

Lengthens the table by 12" (31cm) to 84" (213cm) allowing for access to the patient's neck for carotid artery exams and thyroid exams. 24" (61cm) W x 12" (31cm) L



## QUICKLOCK™ HEAD REST & AERO-CEL PAD #PKG7554T\*


The platform allows precise cervical flexion and height positioning over a wide range for both prone and supine treatments. (Requires Accessory Mounting Plate)



## TERRATOUGH™ FABRIC

Meets the requirements of ISO 10993-1:2009 Biological Evaluation of Medical Tables & passes testing for Cytotoxicity, Skin Irritation & Sensitization as well as CA TB117 Fire Retardance testing. Please refer to the product cleaning instructions for a list of approved cleaners and disinfectants.

## PADDING

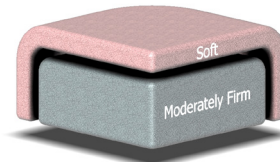

3" Comfort Foam™  
More comfort and support



We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible. However, due to the limitations of printers, scanners, computer monitors, and other equipment, the colors may not be accurate.

CONTACT:

View our full line of  
Medical products

<https://www.oakworksmed.com/pdffdocs/oakworks-medical-catalog/>

717.235.6807 [www.oakworksmed.com](http://www.oakworksmed.com) email: [customerservice@oakworks.com](mailto:customerservice@oakworks.com)