

GENERAL RECTANGULAR TOP Ultrasound Table



Start with a General Ultrasound Table, then simply add the features you want based on the type of ultrasound specialty you perform.

- Flat top table with many optional features to customize your needs
- 550 lb (250kg) bariatric rated lift and load capacity



22" - 38" Height Range



Trendelenburg/Reverse



shown in Sapphire fabric

Inspired by the leading ergonomic experts in sonography



SOUND | ERGONOMICS


Build-A-Bed™ is Trademarked by Sound Ergonomics.

- Decrease sonographer injury
- Features minimize exam times



FM584901

made in the USA
with US & imported parts

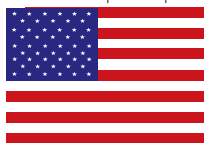

OAKWORKS®
MEDICAL
Position for Success

GENERAL RECTANGULAR TOP Ultrasound Table

Recommended for:

- Abdominal & superficial structures exams
- Musculoskeletal exams
- Vascular exams
- Biopsies



shown in Opal fabric with optional Central Base Locking System and Side Rails

GENERAL RECTANGULAR TOP TABLE

#68666-T* 27" W

#61989-T* 30" W

| FEATURES | |
|-------------------|----------------------------------------------------|
| Rectangular Top | 27" (69cm) or 30" (76cm) wide x 72" (183cm) long |
| Motion 1 | Electric Height Range 22" - 38" (56-97cm) |
| Motion 2 | Electric Trendelenburg/Reverse Tilt $\pm 15^\circ$ |
| Controls | Multi-function hand control |
| Weight Capacity | 550 lb (250kg) |
| Product Weight | 359 lb (163kg) |
| Padding | 3" (8cm) Comfort Foam™ |
| Fabric | Premium, PVC free TerraTouch™ |
| Included Features | 4" (10cm) locking casters, Paper roll holder |
| Warranty | 3 years parts, 2 years labor |

*Fabric color code must be specified.

| OPTIONS must be ordered with table. | | # |
|--------------------------------------------|---------------------------------------------------------------------------|----------|
| Accessory Mounting Plate | Either end | 62123 |
| Adjustable Access Back Support | Left side 0-85° | 73133 |
| Adjustable Access Back Support | Right side 0-85° | 73132 |
| Central Base Locking System | Not available with Patient Lift Compatible Base | 61777 |
| Foot Control | Height, Trendelenburg | 61901 |
| Non-USA electric system | Please contact us for order details | 73134 |
| Patient Lift Compatible Base | Not available with Central Base Locking System | 65893 |
| Side Rails (set of 2) | 23.5" (60cm) long x 6" (15cm) tall (from surface of table) | 61781 |
| Sliding Cardio Panel - Manual | 30" (76cm) width table, left side opens 15" (38cm) | 73130 |
| ACCESSORIES order any time. | | # |
| Detachable Arm Board | Use with T-rail adapter 5" (13cm) wide x 20.5" (52cm) long | 62192 |
| QuickLock™ Head Rest (black) & AeroCel Pad | Comfortable, adjustable prone or supine support (requires mounting plate) | PKG7554* |
| Table Extender | 24" (61cm) wide x 12" (30cm) long (requires mounting plate) | 63487 |

ULTRASOUND FEATURES & OPTIONS

CENTRAL BASE LOCKING SYSTEM #61777

Easily lock your table into place with our unique system that allows locking from both sides of the table.



OPTIONAL

SIDE RAILS #61781

Increase stability and security for your patient. Fold down out of the way when not in use. 23.5" (60cm) long x 6" (15cm) tall (from surface of table)



OPTIONAL

FULL FUNCTION FOOT CONTROLS

#enquire for part number

You can add a foot control that controls all functions.



Multi-function foot control

PATIENT LIFT COMPATIBLE BASE #65893

This option allows the use of patient lifts and easy wheel chair transfers to minimize or eliminate caregiver injuries associated with patient transfers. The high clearance base permits all patient lifts to slide easily under the center of the table for stable transfers while not increasing the low height range of the tables that enable easy wheel chair transfers. Must be ordered with table.

- Includes larger, 5" (13cm) diameter Heavy Duty locking casters
- Provides 6.5" (17cm) clearance for patient lift
- Central Base Locking System not available with this option



OPTIONAL

ADJUSTABLE ACCESS BACK SUPPORT

Left: #73133 Right: #73132

Adjustable with 5 positions 0-85°, allowing easy access for sonographer to the far side and back support for patients. 10" (25cm) wide x 15.5" (39cm) long



OPTIONAL

MANUAL SLIDING CARDIO PANEL

#73130

The sliding cardio panel allows you to adjust the opening to the size you need, up to 15" (38cm). 10" (25cm) wide x 28" (71cm) long. Great for left or right handed scanning. Can also be used as an armrest. Available in manual or electric adjustment options. Available on left side of 30" (76cm) width tables.



OPTIONAL

ACCESSORY MOUNTING PLATE

#62123

Add a face rest or table extender to your table.



OPTIONAL

ULTRASOUND ACCESSORIES

DETACHABLE ARM BOARD #62192-T*

Attaches to the T-rail Adaptor and rotates over 140 degrees. 5" (13cm) wide x 20.5" (52cm) long (Requires T-Rail Adaptor)



T-RAIL ADAPTOR #61902

Multiple locations allow you to choose the best spot for this easily installed & removed accessory capable of handling up to 150 lb (68kg). Can be used with any American standard accessory. 8" (20cm) long

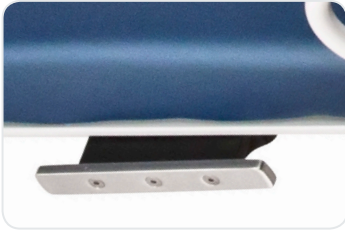


TABLE EXTENDER #63487-T*

Lengthens the table by 12" (31cm) to 84" (213cm) allowing for access to the patient's neck for carotid artery exams and thyroid exams. 24" (61cm) W x 12" (31cm) L



QUICKLOCK™ HEAD REST & AERO-CEL PAD #PKG7554*


The platform allows precise cervical flexion and height positioning over a wide range for both prone and supine treatments. (Requires Accessory Mounting Plate)




TERRATOUGH™ FABRIC

Meets the requirements of ISO 10993-1:2009 Biological Evaluation of Medical Tables & passes testing for Cytotoxicity, Skin Irritation & Sensitization as well as CA TB117 Fire Retardance testing. Please refer to the product cleaning instructions for a list of approved cleaners and disinfectants.

PADDING

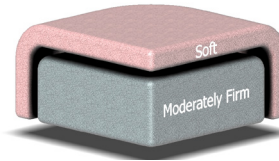

3" Comfort Foam™
More comfort and support

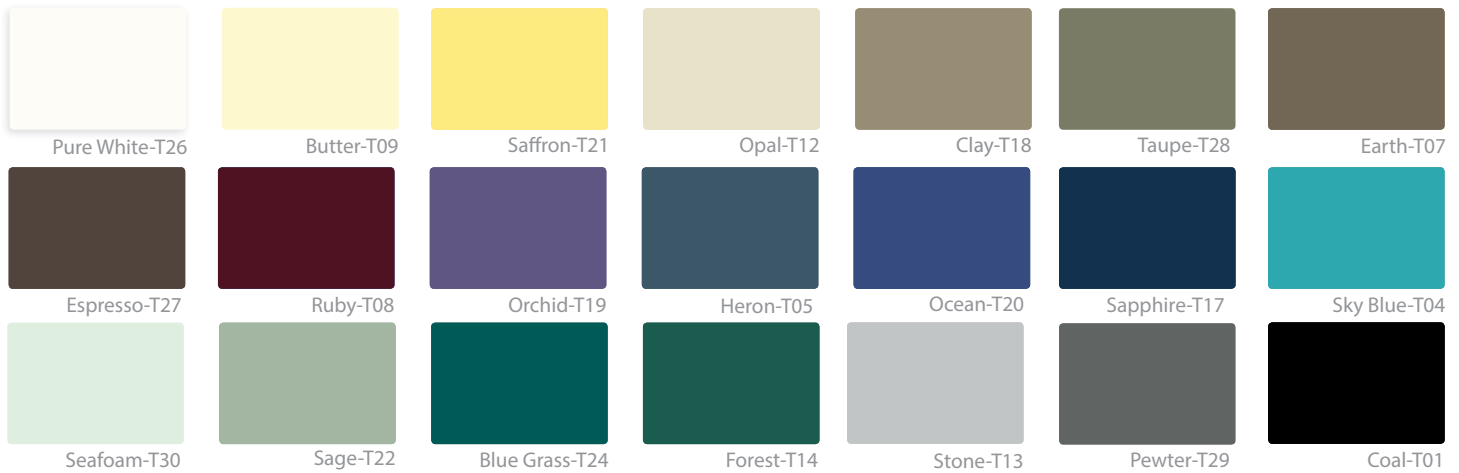

We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible. However, due to the limitations of printers, scanners, computer monitors, and other equipment, the colors may not be accurate.

CONTACT:

View our full line of
Medical products

<https://www.oakworksmed.com/pdfdocs/oakworks-medical-catalog/>

717.235.6807 www.oakworksmed.com email: customerservice@oakworks.com