

# CFPM300 & CFPM301 *Imaging Tables*

- Industry's Largest Adjustable Height Range of 26" - 44"
- 56" Metal Free Imaging Space
- Cost Effective 3 Movement Table
- 500 lb Patient Weight Capacity
- Warranty: 3 years parts, 2 years labor



Integrated Headrest Top



Rectangular Top



HEIGHT RANGE



TRENDELENBURG TILT



LONGITUDINAL TRAVEL  
CFPM301




LATERAL TILT  
CFPM300



FM584901

made in the USA  
with US & imported parts

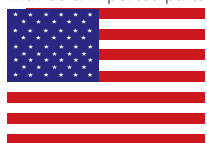

OAKWORKS<sup>®</sup>  
**MEDICAL**  
*Position for Success*

CFPM300 #75208-T01 Integrated  
#75209-T01 Rectangular

CFPM301

#75211-T01 Integrated  
#75212-T01 Rectangular

FEATURES	
Carbon Fiber Top 56" Metal free imaging space	Integrated Headrest 22" (56cm) wide x 84" (214cm) long Rectangular 24" (61cm) wide x 80" (203cm) long
CFPM300 MOVEMENTS	
Motion 1	Electric Height Range 26" - 44" (66-112cm)
Motion 2	Electric Trendelenburg tilt +15°/-12°
Motion 3	Electric Lateral tilt ±15°
CFPM301 MOVEMENTS	
Motion 1	Electric Height Range 26" - 44" (66-112cm)
Motion 2	Electric Trendelenburg tilt +15°/-11°
Motion 3	Electric Longitudinal travel 10" (25cm)
Controls	Multi-function hand control
Weight Capacity	500 lb (227kg)
Product Weight	475 lb (216kg)
Included Features	1" (3cm) Radiolucent Comfort Foam™ table top pad, 4" (10cm) locking casters, Crescent face rest pad (included with integrated headrest top), patient safety straps (2)
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, CE marked, cETLus listed, CB Certificate

OPTIONS <i>must be ordered with table.</i>		#
2" (5cm) Table Top Pad Upgrade	Upgrade from the standard 1" (3cm) pad	73190-T01
5" (13cm) casters	Upgrade from the standard 4" (10cm)	65427
Foot Control (CFPM300)	multi-functional with cradle on table base	75074
Foot Control (CFPM301)	multi-functional with cradle on table base	75076
Non-USA electric system	Please contact us for order details	73134
ACCESSORIES <i>order any time. See accessories on pages 34-39.</i>		#
Arm Hammock	Use with Integrated Head Rest Top	48751-T01
Carbon Fiber Arm Board (6" wide)	180° Articulation	42585
Carbon Fiber Arm Board (7.5" wide)	180° Articulation	70337
Fluoro Extender	Polycarbonate	33975
Interventional Supine Pillow System	Contoured pad and support wedge pad	69808-T01
Spine Positioning System II	8 piece system	70100-T01
T-rail Clamp	6" (15cm) long rail	63336

## SINGLE CARBON FIBER ARM BOARD

#42585-T01 - 6" wide #70337-T01 - 7.5" wide

Fully radiolucent with a 180° range of motion, the arm board requires no mounting or attachment hardware and simply slides under the table top pad. It can be used on any fluoroscopy table and will support arms and feet in a variety of positioning configurations. Folds for storage.



## SPINE POSITIONING SYSTEM II #70100-T01

- **Optimal Imaging**  
High radiolucency and flexibility for ideal cervical, thoracic and lumbar visualization.
- **Enhanced Patient Comfort**  
Minimize unwanted movement and reduce procedure time.
- **Stable Design**  
Large contoured torso pad cradles patient & the strap securely fastens the system to the table.



## 6" ADJUSTABLE T-RAIL CLAMP #63336

The adjustable T-Rail clamp can be positioned along the grooved table top edge. Clamp is also removable.

- Accommodates a variety of T-rail attachable accessories
- 150 lb (68kg) weight rating



CONTACT:

View our full line of  
Medical products



Scan this code with  
your smart phone

#72703



717.235.6807 www.oakworksmed.com email: customerservice@oakworks.com